

2.1.3 Install USB3.0 driver

When using USB 3.0 download with ACS camera, USB 3.0 driver must be installed.

Open the folder corresponding to your Windows from the "USB 3.0 Drivers" folder on the installation CD. (Example: Windows 10 64bit environment: USB 3.0 Drivers ¥ Win10 ¥ x64)

Please execute the installer in the folder, "DPinst32.exe" for 32bit environment,

"DPinst64.exe" for 64bit environment. If you are asked to confirm user account control, click "Yes" and install according to the instructions of the wizard.

Click the "Finish" button to close the wizard.

2.1.4 Install .NET Framework 4.7.1

“.NET Framework version 4.7.1” or later must be installed for operating all function of MLink.


If “.NET Framework version 4.7.1” or later is not installed on the Windows used, execute

“NDP471-KB4033342-x86-x64-AllOS-ENU.exe” in the “MLink” folder on the installation CD to install.

To install ".NET Framework 4.7.1", the OS must be Windows 7 SP1 or later version.

2.1.5 Setting IP Address on Installing PC

Open the control panel from the Windows start button.

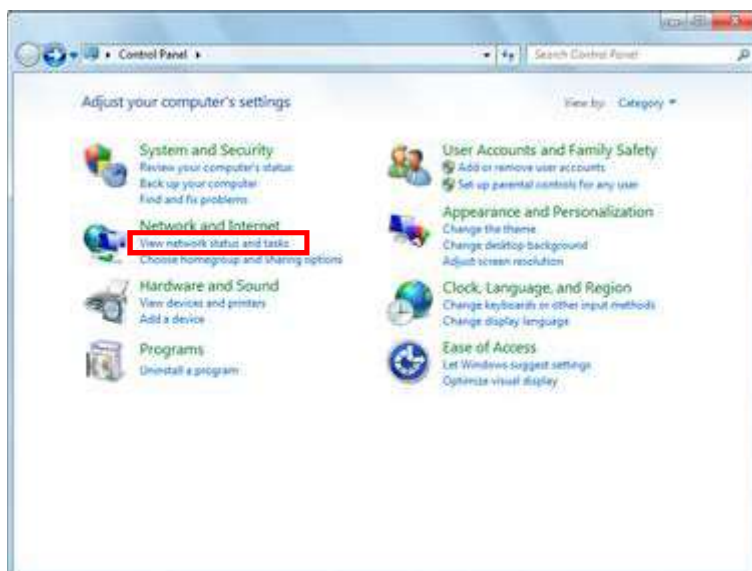

Click View network status and tasks.